

PRODUCT DATA SHEET

Bioworld Technology,Inc.

mGLUR2 polyclonal antibody

Catalog: BS67494 Host: Rabbit Reactivity: Human, Mouse, Rat

BackGround:

Metabotropic glutamate receptor 2 (mGluR2) is a class C G protein-coupled receptor for the neurotransmitter glutamate in the mammalian brain. Unlike ionotropic receptors, metabotropic receptors do not form an ion channel pore themselves but are indirectly linked to ion channels . While mGluR1 and mGluR5 activate phospholipase C, mGluR2, mGluR3, mGluR4, and mGluR6 are coupled to the inhibitory G protein $G\alpha(i/o)$ and inhibit adenylyl cyclase (AC) activity . Research studies have suggested that mGluR2/3 receptors may be potential targets for the treatment of Schizophrenia . Furthermore, mGluR2 interacts with the 5HT2A serotonin receptor to form a hetero-complex in the brain. This complex is a potential pharmacological target for hallucinogenic drugs .

Product:

Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol, and 0.01% sodium azide.

Molecular Weight:

~ 100 kDa

Swiss-Prot:

Q14416

Purification&Purity:

The antibody was purified by immunogen affinity chromatography.

Applications:

WB (1/500 - 1/1000), IHC (1/50 - 1/200), IF/ICC (1/50 - 1/200)


Storage&Stability:

Store at 4 ${\mathbb C}$ short term. Aliquot and store at -20 ${\mathbb C}$ long term. Avoid freeze-thaw cycles.


Specificity:

Recognizes endogenous levels of mGLUR2 protein.


DATA:

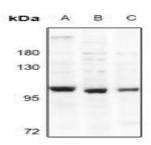

Western blot analysis of mGLUR2 expression in U87MG (A), A375 (B), C6 (C) whole cell lysates.

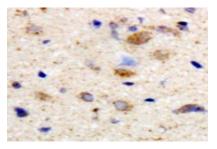

Immunohistochemical analysis of mGLUR2 staining in human brain formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

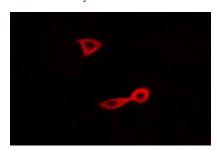

Immunofluorescent analysis of mGLUR2 staining in LOVO cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a hidified chamber. Cells were washed with PBST and incubated with a Alexa Fluor

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841

Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151



PRODUCT DATA SHEET

Bioworld Technology,Inc.

594-conjugated secondary antibody (red) in PBS at room temperature in the dark.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>
Tel: 6123263284

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151