

## PRODUCT DATA SHEET

Bioworld Technology,Inc.

# Peroxiredoxin 1 monoclonal antibody

Catalog: MB65970 Host: Mouse Reactivity: Human, Mouse, Rat

#### **BackGround:**

Thiol-specific peroxidase that catalyzes the reduction of hydrogen peroxide and organic hydroperoxides to water and alcohols, respectively. Plays a role in cell protection against oxidative stress by detoxifying peroxides and as sensor of hydrogen peroxide-mediated signaling events. Might participate in the signaling cascades of growth factors and tumor necrosis factor-alpha by regulating the intracellular concentrations of H2O2 (PubMed:9497357).

Reduces an intramolecular disulfide bond in GDPD5 that gates the ability to GDPD5 to drive postmitotic motor neuron differentiation (By similarity).

#### **Product:**

Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol, and 0.01% sodium azide.

#### **Molecular Weight:**

~ 21 kDa

## **Swiss-Prot:**

Q06830

## **Purification&Purity:**

The antibody was affinity-purified from mouse antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

#### **Applications:**

WB (1/1000 - 1/3000), IF/ICC (1/100 - 1/200)


#### **Storage&Stability:**

Store at  $4\,\mathrm{C}$  short term. Aliquot and store at  $-20\,\mathrm{C}$  long term. Avoid freeze-thaw cycles.


## **Specificity:**

Recognizes endogenous levels of Peroxiredoxin 1 protein.

## **DATA:**

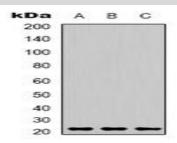

Western blot analysis of Peroxiredoxin 1 expression in MCF7 (A), mouse brain (B), rat kidney (C) whole cell lysates.

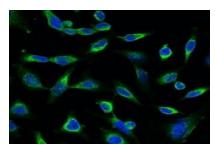

Immunofluorescent analysis of Peroxiredoxin 1 staining in Hela cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a hidified chamber. Cells were washed with PBST and incubated with a FITC-conjugated secondary antibody (green) in PBS at room temperature in the dark. DAPI was used to stain the cell nuclei (blue).

#### Note:

For research use only, not for use in diagnostic procedure.

## Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151